

# CREATININE-E FL

CE F125 CH	5 x 25 ml
CE F375 CH	15 x 25 ml
CE F600 CH	10 x 60 ml

## INTENDED USE

Reagent for quantitative *in vitro* determination of creatinine in biological fluids.

## SUMMARY OF TEST

Between 1 and 2% of muscle creatine is converted to creatinine daily. Because the amount of endogenous creatinine produced is proportional to muscle mass, the production varies with age and sex. Because creatinine is endogenously produced and released into body fluids at a constant rate and its plasma levels maintained within narrow limits, its clearance may be measured as an indicator of glomerular filtration rate (GFR).

## PRINCIPLE OF THE METHOD

Through a series of enzymatic reactions, creatinine is converted in glycine, whilst endogenous components such as creatine and sarcosine are eliminated in the first step of the sequence. The formed hydrogen peroxide reacts with TOPS in the presence of peroxidase, to give a quinoneimine dye. The intensity of color, measured at 546 nm, is proportional to creatinine concentration in the sample.

## KIT COMPONENTS

### For *in vitro* diagnostic use only.

The components of the kit are stable until expiration date on the label at 2-8°C.

Keep away from direct light sources.

**CREA-E R1** F125: 4 x 25 ml (liquid) blue cap  
F375: 12 x 25 ml (liquid) blue cap  
F600: 8 x 60 ml (liquid) blue cap

**CREA-E R2** F125: 1 x 25 ml (liquid) red cap  
F375: 3 x 25 ml (liquid) red cap  
F600: 2 x 60 ml (liquid) red cap

Composition in the test: Creatinase  $\geq$  10 kU/l, Creatininase  $\geq$  10 kU/l, Sarcosine Oxidase  $\geq$  1 kU/l, Peroxidase  $\geq$  5 kU/l, TOPS  $\geq$  3 mM, 4-aminoantipyrine  $\geq$  20 mg/l.

Store all components at 2-8°C.

## MATERIALS REQUIRED BUT NOT SUPPLIED

Current laboratory instrumentation. Spectrophotometer UV/VIS with thermostatic cuvette holder. Automatic micropipettes. Glass or high quality polystyrene cuvettes. Saline solution.

## REAGENT PREPARATION

Use separate reagents ready to use.

Stability: up to expiration date on labels at 2-8°C.

Stability since first opening of vials: use preferably within 60 days at 2-8°C

## PRECAUTIONS

Reagent may contain some non-reactive and preservative components. It is suggested to handle carefully it, avoiding contact with skin and swallow.

Perform the test according to the general "Good Laboratory Practice" (GLP) guidelines.

## SPECIMEN

Serum, plasma. Urine.

Creatinine is stable 24 hours at 2-8°C. Freeze samples for prolonged storage.

Dilute urine sample 1:100 with deionized water.

## TEST PROCEDURE

Wavelength:	546 nm		
Lightpath:	1 cm		
Temperature:	37°C		
dispense:	blank	calibrator	sample
reagent R1	1 ml	1 ml	1 ml
water	20 $\mu$ l	-	-
calibrator	-	20 $\mu$ l	-
sample	-	-	20 $\mu$ l
Mix, incubate at 37°C for 5 minutes. Read absorbances of calibrator ( $Ac_1$ ) and samples ( $Ax_1$ ) against reagent blank.			
dispense:	blank	calibrator	sample
reagent R2	250 $\mu$ l	250 $\mu$ l	250 $\mu$ l
Mix, incubate at 37°C for 5 minutes. Read absorbances of calibrator ( $Ac_2$ ) and samples ( $Ax_2$ ) against reagent blank.			

## RESULTS CALCULATION

Serum/plasma sample:

$$\text{creatinine mg/dl} = (Ax_2 - Ax_1) / (Ac_2 - Ac_1) \times \text{calibrator value}$$

Random urine sample:

$$\text{creatinine mg/dl} = (Ax_2 - Ax_1) / (Ac_2 - Ac_1) \times \text{calibrator value} \times 10 \text{ (dilution)}$$

24 hours urine sample (creatinine mg/24h):

$$\text{creatinine mg/24h} = (Ax_2 - Ax_1) / (Ac_2 - Ac_1) \times \text{calibrator value} \times 10 \times \text{diuresis (dilution, diuresis in dl)}$$

## EXPECTED VALUES

Serum/plasma:

Men: 0.67 - 1.17 mg/dl (59 - 104  $\mu$ mol/l)

Women: 0.51 - 0.95 mg/dl (45 - 84  $\mu$ mol/l)

24h urine:

Men: 1000 - 2000 mg/24h (8.85 - 17.70 mmol/24h)

Women: 800 - 1800 mg/24h (7.08 - 15.93 mmol/24h)

Each laboratory should establish appropriate reference intervals related to its population.

## QUALITY CONTROL AND CALIBRATION

It is suggested to perform an internal quality control. For this purpose the following human based control sera are available:

### QUANTINORM CHEMA

with normal or close to normal control values.

### QUANTIPATH CHEMA

with pathological control values.

If required, a multiparametric, human based calibrator is available:

### AUTOCAL H

Please contact Customer Care for further information.

## TEST PERFORMANCE

### Linearity

The method is linear up to 50 mg/dl.

If the value is exceeded, it is suggested to dilute sample 1+9 with saline and to repeat the test, multiplying the result by 10.

### Sensitivity/limit of detection (LOD)

The limit of detection is 0.04 mg/dl.

### Interferences

No interference was observed by the presence of:

hemoglobin  $\leq$  1000 mg/dl

bilirubin  $\leq$  28 mg/dl

lipids  $\leq$  1400 mg/dl

ascorbic acid  $\leq$  50 mg/dl

### Precision

intra-assay (n=10)	mean (mg/dl)	SD (mg/dl)	CV%
sample 1	0.98	0.01	1.10
sample 2	3.98	0.02	0.56

inter-assay (n=20)	mean (mg/dl)	SD (mg/dl)	CV%
sample 1	0.97	0.03	3.28
sample 2	3.97	0.11	2.90

## Methods comparison

A comparison between Chema and a commercially available product gave the following results with 99 samples:

$$\text{Creatinine competitor} = x$$

$$\text{Creatinine Chema} = y$$

$$y = 1.004x + 0.037 \text{ mg/dl} \quad r^2 = 0.998$$

## WASTE DISPOSAL

This product is made to be used in professional laboratories.

P501: Dispose of contents according to national/international regulations.







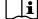
## REFERENCES

Tietz Textbook of Clinical Chemistry, Fourth Edition, Burtis-Ashwood-Bruns (2006), 797-801  
Clin. Chem. 2012, 58(2), 391-401

## MANUFACTURER

Chema Diagnostica  
Via Campania 2/4  
60030 Monsano (AN) - ITALY - EU  
phone +39 0731 605064  
fax +39 0731 605672  
e-mail: mail@chema.com  
website: http://www.chema.com

## SYMBOLS

	<i>in vitro</i> diagnostic medical device
	batch code
	catalogue number
	temperature limit
	use by date
	caution
	consult instructions for use